### **OptiView**. Access. Comfort. Even retraction.



# LIP AND CHEEK RETRACTOR













#### NEW! SMALL SIZE

## Enjoy the unobstructed view.

- Superior, uniform perioral retraction. Maximizes visibility and access to the working field.
- Passive placement. Stable, flexible 3-dimensional design ensures quick and gentle placement and removal.
- Unobstructed access. No additional hands or instruments required.
- Enhanced patient comfort. Anatomic design ensures wearing comfort even during longer procedures and facilitates placement of the saliva ejector.
- Dual-duty. Provides both lip and cheek retraction with one device.




#### OptiView. Total access and visibility with no extra hands.

OptiView is an innovative device that provides even perioral retraction of the lips and cheeks to enlarge the treatment field for greatly improved access, visibility and patient comfort.

Most conventional retractors provide either lip or cheek retraction, but not both. And due to their high structural tension, they tend to be uncomfortable for the patient. Since OptiView conforms to the anatomy of the mouth and minimizes tension, it is comfortable for patients to wear for longer periods than was ever possible.

**Indications:** Suitable for use in diagnostic, preventive and therapeutic dental applications such as intraoral photography, tooth preparation, impression taking, oral prophylaxis, bleaching and occlusal analysis.

# Optimal perioral retraction without requiring instruments or a second pair of hands to hold the device in place.

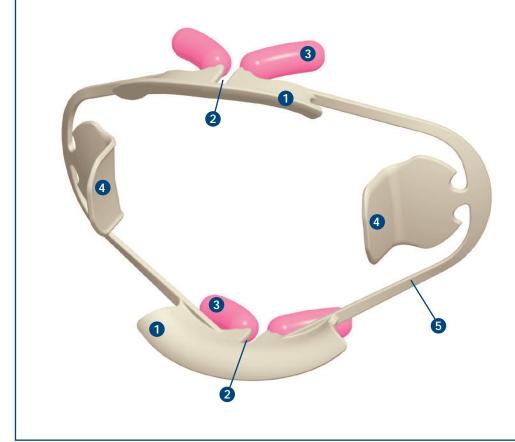

- Anatomically shaped maxillary and mandibular lip holders provide superior access to the buccal and gingival areas.
- 2. Lip holders include opening for the frenulum to maximize comfort.
- 3. Cushioned lip supports position lips an ample distance away from anterior teeth and gums.
  - Autoclavable cushions can be removed from the frame to facilitate cleaning.
- Lateral flaps gently hold cheeks and contract for easy placement and removal.
- High-tech OptiView material is safe to autoclave at 134°C/273°F (3 minutes).



#### **Order Information**

5500 - OptiView Standard

Contains: 2 Lip and Cheek Retractors with Mounted Cushions

5502 - OptiView Small

Contains: 2 Lip and Cheek Retractors Small size with Mounted Cushions

5501 - OptiView Refill Cushions (pack of 10)

